
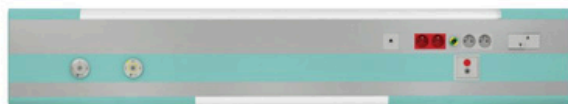


MERYlight is a very functional bed head unit designed for modern hospital rooms. It provides convenient and safe access to all necessary fittings, combining gas, electrical & IT sockets, as well as lighting systems in one compact design.

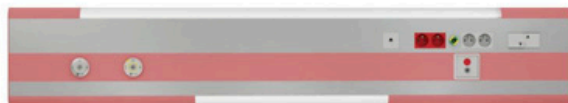

TECHNICAL SPECIFICATION

- Body made of anodized aluminium profile, mounted directly on the wall
- Gas front cover powder coated in any colour (RAL colour palette)
- Number of inlets (G/ E) depending on configuration and fittings

Size:

- Length (standard): single unit-1600 mm, double unit-3200 mm, triple unit-4800 mm
- Non-standard lengths available
- Depth: 65 mm
- Height: 290 / 300 / 305 / 310 mm – depending on configuration

MEDICAL GAS OUTLETS AND VACUUM

- Available standards: SS / DIN / BS / AFNOR / other
- Medical gases location on the BHU: at the door side

LIGHTING

- Types of lighting: LED / 2G11
- Lighting configuration: general/room, reading, night
- Lighting is switched on/off by switch/button on the BHU and/or the wall and/or patient's handset

CHARACTERISTICS

- The slimmest BHU in INMED range – depth 65mm
- Medical gas outlets, electrical and IT sockets fitted in front panel



STANDARD CONFIGURATION IN HOSPITAL WARDS

- Single unit 1600 mm: recommended
- Medical gas outlets: 1x O₂, 1x VAC
- Electrical sockets: 230 V, 3 pcs (2 circuits)
- Equipotential pins: 1 pc
- IT sockets: RJ-45, 6 cat., 1 pc
- General lighting 1x 24W LED, 4400 lm: switch on the BHU
- Reading lighting 1x 12W LED, 2200 lm: switch on the BHU
- Night lighting 1x 4W LED, 300 lm: switch on the BHU
- NCS cut-out

ELECTRICAL AND IT SOCKETS

- Electrical sockets: Type E (FR) / Type F (Schuko) / Type G (BS) / other
- IT sockets: RJ-11 / RJ-12 / RJ-45, 6a cat. or higher
- Multimedia sockets: USB, HDMI, and other
- NCS depending on customers' demand

VERSION A

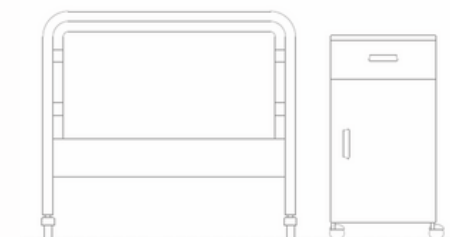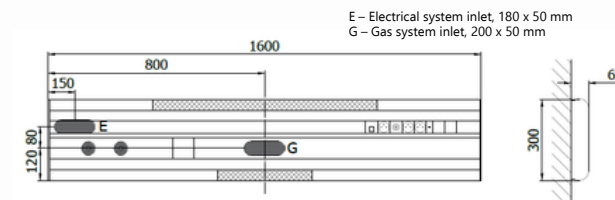

VERSION B

(mirror reflection)

