

MERYart is a line of modern premium bed head units where medical functionality and high design esthetics come together. Available with Corian® or tempered glass fronts, MERYart accommodates medical gases, electrical and IT systems. Customisable graphics are available.



TECHNICAL SPECIFICATION

- Body made of anodized or powder coated aluminium profile mounted directly on the wall
- Decorative front panel made of tempered glass

Size:

- Length (standard): single unit – 1500 mm,
- Depth: 124 mm
- Height: 270 mm

MEDICAL GAS OUTLETS AND VACUUM

- Available standards: SS / DIN / BS / AFNOR / other
- Medical gases location on the BHU: at the door side

LIGHTING

- Types of lighting: LED / 2G11
- Lighting configuration: general/room, reading, night
- Lighting is switched on/off by switch/button on the BHU and/or the wall and/or patient's handset

CHARACTERISTICS

- Electrical and IT sockets fixed at an angle of 25° to the floor as in MERYclassic type
- Graphics selection customized depending on customers demand



STANDARD CONFIGURATION IN HOSPITAL WARDS

- Medical gas outlets: 1x O2, 1x VAC
- Electrical sockets 230 V: 3 pcs (2 circuits)
- Equipotential pins: 1 pc
- IT sockets: RJ-45, cat. 6: 1 pc
- Reading lighting 1x 12W LED, 2200 lm: switch on the unit
- Night lighting 1x 4W LED, 300 lm: switch on the unit
- NCS cut-out

ELECTRICAL AND IT SOCKETS

- Electrical sockets: Type E (FR) / Type F (Schuko) / Type G (BS) / other
- IT sockets: RJ-11 / RJ-12 / RJ-45, 6a cat. or higher
- Multimedia sockets: USB, HDMI, other
- NCS depending on customers' demand

VERSION A

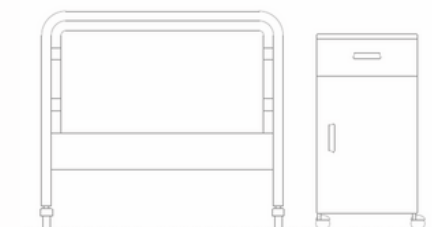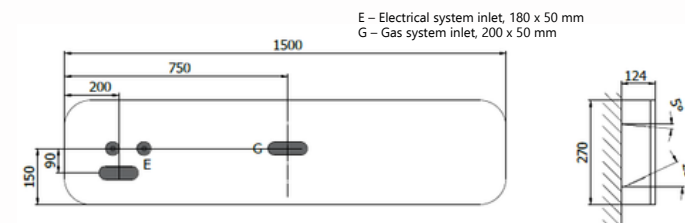

VERSION B

(mirror reflection)

