


MERYedge is a modern bed head unit coming with elegant brushed aluminium front, designed for functionality and esthetics of hospital rooms. It is a unit where all media, advanced lighting and well thought-out ergonomics come together.

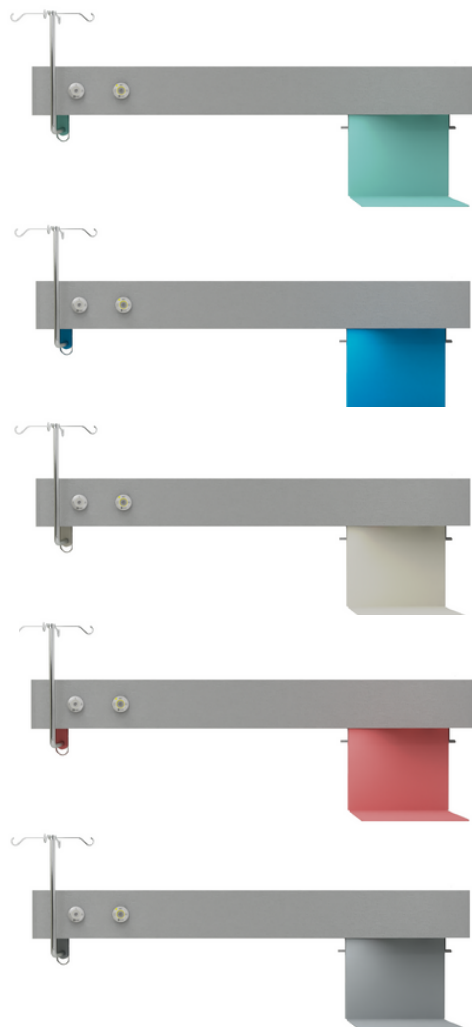

TECHNICAL SPECIFICATION

- Body made of anodized or powder-coated aluminium profile mounted directly on the wall
 - Front cover made of anodized aluminum, brushed
 - Shelf and rail fixed on the unit as option
- Size:
- Length (standard): single unit – 1500 mm, double unit – 3200 mm, triple unit – 4800 mm
 - Depth: 160 mm
 - Height: 164 mm

MEDICAL GAS OUTLETS AND VACUUM

- Available standards: SS / DIN / BS / AFNOR / other
- Medical gases location on the BHU: at the door side

LIGHTING

- Lighting types: LED / 2G11
- Lighting configuration: general/room, reading, night
- Lighting is switched on/off by switch/button on the BHU and/or the wall and/or patient's handset

CHARACTERISTICS

- Electrical sockets and lighting switches fitted in the bottom of the unit, parallel to the ground
- Front panel coming as one block with end sides

STANDARD CONFIGURATION IN HOSPITAL WARDS



- Medical gas outlets: 1x O2, 1x VAC
- 230 V electrical sockets: 3 pcs (2 circuits)
- Equipotential pins: 1 pc
- IT sockets RJ-45, cat. 6: 1 pc
- General lighting 1x 24W LED, 4400 lm: switch on the unit
- Reading lighting 1x 12W LED, 2200 lm: switch on the unit
- Night lighting 1x 4W LED, 300 lm: switch on the unit
- NCS cut-out

ELECTRICAL AND IT SOCKETS

- Electrical sockets: Type E (FR) / Type F (Schuko) / Type G (BS) / other
- IT sockets: RJ-11 / RJ-12 / RJ-45, 6a cat. or higher
- Multimedia sockets: USB, HDMI, other
- NCS depending on customers' demand

VERSION A

E – Electrical system inlet, 180 x 50 mm
G – Gas system inlet, 200 x 50 mm

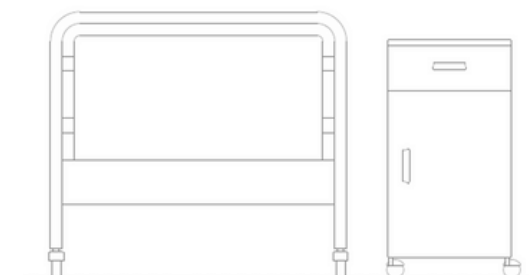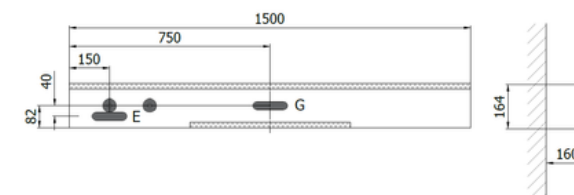

VERSION B

(mirror reflection)

