



**LISSAclassic-lift** ceiling pendant with lifting arm is designed to gas and electricity distribution and provides ergonomic access to equipment necessary for the work of medical personnel. Due to service head lifting, **LISSAclassic-lift** pendant enables precise adjustment of working height to the user's needs, ensuring full ergonomics and safety at work.


This pendant is used in:

- operating rooms
- resuscitation rooms
- treatment rooms and other rooms where workspace optimisation is key.

Depending on configuration LISSAclassic-lift can function as (among others):

- **anaesthesia pendant**
- **surgical pendant**
- **auxiliary pendant in operating rooms.**

It is particularly recommended for use as anaesthesia pendant, used with mobile anaesthesia machine, enabling convenient and safe management of all necessary connections without restricting equipment mobility.

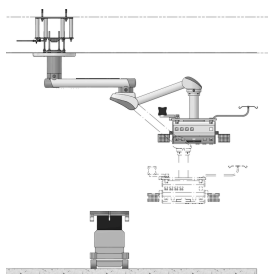

### ARMS TECHNICAL SPECIFICATION

- Light arm
- Single/double swivel arm
- Single arm length: 1000mm
- Total arms range: 1600, 1800, 2000, 2200, 2400 mm
- Possibility of using non-standard arms ranges
- Rotation – maximum rotation range: 350°, can be reduced by adjusting rotation limiters
- Brake: pneumatic, electromagnetic or friction
- Arms, ceiling covers, service head covers and accessories powder coated in RAL 9002
- Arms end-caps, drop tubes, bottom covers and blind covers powder coated in RAL 9022

### SERVICE HEAD TECHNICAL SPECIFICATION

- Horizontal service heads width: 600, 700, 800, 900, 1000 mm
- Non-standard size horizontal and vertical service heads possible
- Service head's body made of aluminium profiles
- Non-standard RAL for service head covers possible

### MEDICAL GAS AND VACUUM OUTLETS

- Standards: SS / DIN / BS / AFNOR / other
  - Type: O<sub>2</sub>, AIR, N<sub>2</sub>O, CO<sub>2</sub>, VAC
- Gas outlets quantities upon to individual configuration
- AGSS outlet

### ELECTRICAL AND IT SOCKETS

- Electrical sockets: Type E (FR) / Type F (Schuko) / Type G (BS) / other
- IT sockets: RJ-11 / RJ-12 / RJ-45, 6a cat. or higher
- Multimedia sockets: USB, HDMI and other
- Nurse call system upon to user's needs

### OPTIONAL LIGHTING

- in upper part of ceiling arm
- in side cover of service head
- in the bottom of service head
- mobile (individual lamp), fixed on the rail

### OPTIONAL EQUIPMENT

- Medical rails made of acid-resistant steel
- Additional equipment:
  - IV poles
  - monitors and patient monitors hangers / arms
  - medical disposables baskets and other

#### Vertical head sections

Rectangular Square Trapezoidal

