



VEGAcare ceiling-floor medical supply unit is intended to gas and electricity distribution and provides ergonomic access to equipment necessary for the work of medical personnel. It has been designed for use in rooms where ceiling mounting is impossible, difficult or not recommended for structural or installation reasons.


Thanks to its brace load-bearing construction, **VEGAcare** is assembled directly between floor and ceiling, eliminating the need to interfere with building structure or reinforce the ceiling. VEGAcare provides all functionality typical for ceiling pendants, while maintaining ease of installation and exceptional stability.

VEGAcare is used in:

- **Intensive Care Units**
- **Treatment and dressing rooms**
- **Hospital rooms with lighter equipment requirements**
- **Auxiliary and preparation rooms**
- **Hospital emergency departments**

Depending on configuration VEGAcare can function as:

- **ICU pendant**
- **anaesthesia pendant**
- **treatment pendant**

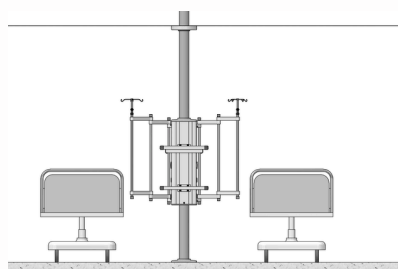

SERVICE HEAD TECHNICAL SPECIFICATION

- Vertical 3, 4 or 6-plane service heads length: 600, 800, 1000, 1200, 1500 mm
- Non-standard size service heads possible
- Service head's body made of aluminium profiles
- Ceiling covers, service head covers and accessories powder coated in RAL 9002
- Bearings end-caps, drop tubes, bottom covers and blind covers powder coated in RAL 9022
- Non-standard RAL for service head covers possible

MEDICAL GAS AND VACUUM OUTLETS

- Standards: SS / DIN / BS / AFNOR / other
 - Type: O₂, AIR, N₂O, CO₂, VAC
- Gas outlets quantities upon to individual configuration
- AGSS outlet

ELECTRICAL AND IT SOCKETS

- Electrical sockets: Type E (FR) / Type F (Schuko) / Type G (BS) / other
- IT sockets: RJ-11 / RJ-12 / RJ-45, 6a cat. or higher
- Multimedia sockets: USB, HDMI and other
- Nurse call system upon to user's needs

OPTIONAL LIGHTING

- in side cover of service head
- in the bottom of service head
- mobile (individual lamp), fixed on the rail

OPTIONAL EQUIPMENT

- Vertical poles for accessories
- Medical rails made of acid-resistant steel
- Different size and load capacity shelves
- Drawers
- IV poles
- Monitor hangers / arms
- Catheters and medical disposables baskets
- Cables and hoses organizers

Vertical head sections

